




AUTOPSY TABLE

Autopsy Table

Autopsy Table is the series of Mortuary and histology products. Unit is fabricated from the high quality gauge stainless steel 304. It is with the fashion and the humanized design. For using in histology lab, Human Anatomy lab, and Forensic appraisal department. In order to meet the different requirements of users, we support several models; Simple model and the luxury models. Safe and high quality work will bring users the effective.



SPINNE1

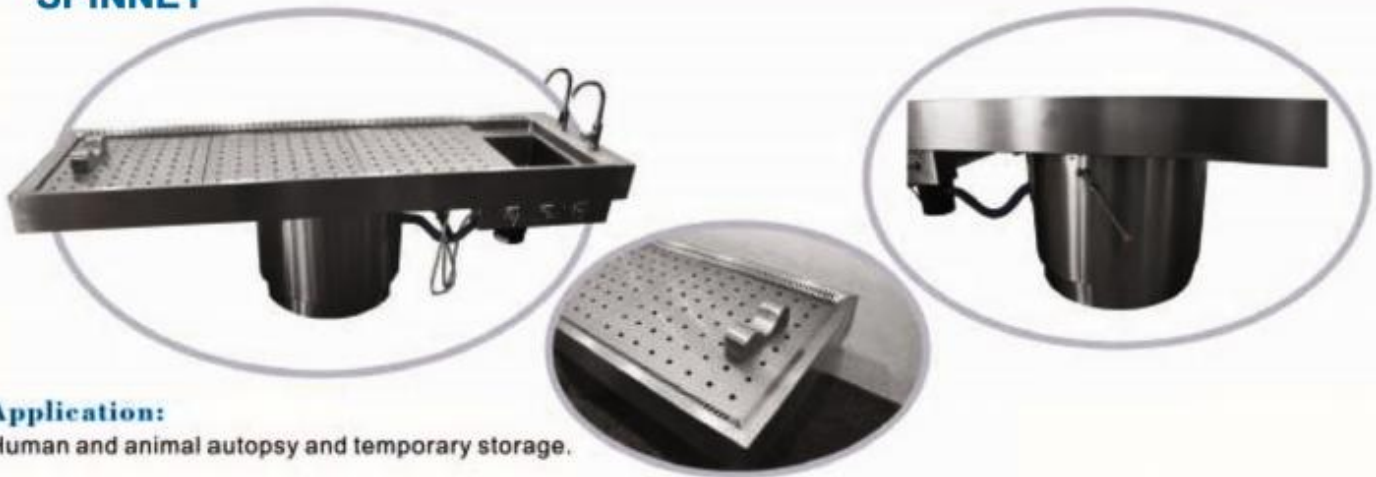

Application:

Human and animal autopsy and temporary storage.

Feature

All stainless steel 304 construction. It is half-closed structure.

Heliarc-welded seams and joints, ground and polishes to match adjoining surfaces.

It with the adjustable square base. The table height is adjustable from 720 to 920mm.

Table can rotate 40 degrees.

Operating tray Removable(3pcs) and detachable operating tray.

Plumbing and electrical lines factory installed and require only single point connection at the project location.

There is the the power of the socket on the side of table. It is waterproof . Can be used for other instrument connections.

Crushing and sewage system:

Large single compartment sink basin.

With the pulverizer, It can crush dirt and little bone instantly.

Two taps. Sensor faucet with cold water and manual faucet with hot or cold water.

With the removable spray tap, the water pressure can be adjusted. It can wash and clean each corner of table dirt

Drained system

Rapid and positive drainage

With drained pipeline, the treated liquid drained along the drained system.



Exhaust system function. Air outlet is at the inside of the table. The exhaust way is down. It needs to connect with the exhaust blower outside.

OEM and customer special request is accepted.

Specification

Model Number	SPINNE1
Material	Stainless Steel 304
Product Size (mm)	L2600 x W800 x H 800
	L102.3 x W31.5 x H 31.5
Table height(adjustable)	720—920mm (Optional); By motor.
Table can rotate	40 degrees; Manual.
Voltage	AC 220v/50hz; electricity 1.5A;
Power for the Pulverizer	220v/50hz; 370w; frequency : 1700times/minute;
Load Capacity	200KG
Package Size	2800*900*900
Table height(adjustable)	110.2*35.4*35.4
Volume(m3)	2.27CBM
Package	Standard exporting, plywood case.

SPINNE2

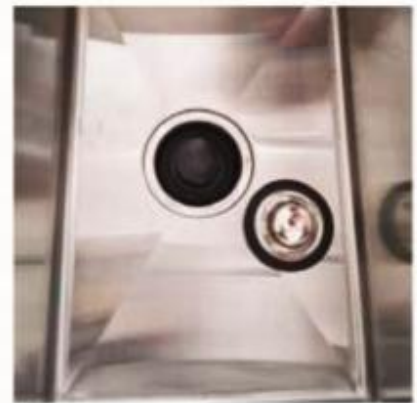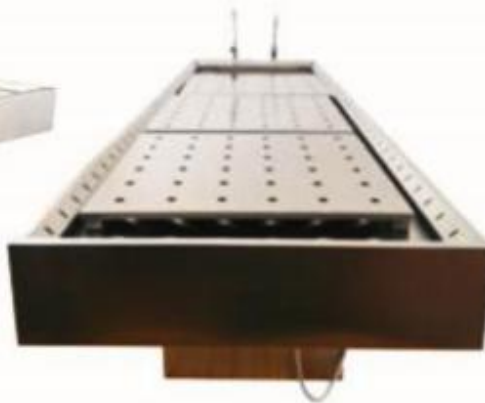

Application:

Human and animal autopsy and temporary storage.

Feature

All stainless steel 304 construction. It is half-closed structure.

Heliarc-welded seams and joints, ground and polishes to match adjoining surfaces.

It with the adjustable square base. The table height is adjustable from 720 to 920mm (Optional) .

Operating tray Removable(3pcs) and detachable operating tray.

Plumbing and electrical lines factory installed and require only single point connection at the project location.

There is the the power of the socket on the side of table. It is waterproof . Can be used for other instrument connections.

Crushing and sewage system:

Large single compartment sink basin.

With the pulverizer, It can crush dirt and little bone instantly.

Two taps. Sensor faucet with cold water and manual faucet with hot or cold water.

With the removable spray tap, the water pressure can be adjusted. It can wash and clean each corner of table dirt

Drained system

Rapid and positive drainage

With drained pipeline, the treated liquid drained along the drained system.

Exhaust system function. Air outlet is at the inside of the table, The exhaust way is down. It needs to connect with the exhaust blower outside.

OEM and customer special request is accepted.

Specification

Model Number	SPINNE2 Fix height model
	SPINNE2 Height adjustable model
Material	Stainless Steel 304
Product Size (mm)	L2600 x W800 x H 800
	L102.3 x W31.5 x H 31.5
Table height(adjustable)	720—920mm (Optional); By motor. Optional
Voltage	AC 220v/50hz; electricity 1.5A;
Power for the Pulverizer	220v/50hz; 370w; frequency : 1700times/minute;
Load Capacity	200KG
Package Size	2800*900*900
	110.2*35.4*35.4
Volume(m ³)	2.27CBM
Package	Standard exporting, plywood case.

SPINNE3



Application

Human and animal autopsy and temporary storage

Feature

All stainless steel 304 construction.

Heliarc-welded seams and joints, ground and polishes to match adjoining surfaces

There is a hole at the middle of table, to be easy for cleaning table.

OEM and customer special request is accepted.

Drained system

With drained pipeline, the treated liquid drained along the drained system.

Removable shower head.

The pipe of removable shower head is long enough to clean the whole operating table body.

Package information

Package Standard export package: plywood case

Specification

Model Number	SPINNE3
Material	Full stainless steel 304
Cleaning system	Cleaning system
	Removable tap
Product Size (mm)	L2300*W800*H800
Table Height (mm)	800
Water bath size (mm)	L380*W320*H150
Package Size (mm)	2400*900*600
Volume(m ³)	1.20
Load Capacity	200kg
Package	Standard exporting, plywood case.

SPINNE4

Product Name:

Anatomical / autopsy table

Applications:

Mortuary; Anatomy, forensic examinations, medical colleges and hospitals, etc.

Product size:

L2500*W800*H750 to 900mm; Special sizes can be customized.

Package:

Packed in plywood case; packing size:2630*850*880MM

Specification:

- Overall 304 stainless steel structure. Surface grinding and wire drawing.
- All internal corner arc design for easy cleaning.
- The table can be manually rotated by 45° and it can be fixed at any angle.
- 3 perforated, removable panels.
- Ventilation: Downward suction.
- Aspirator and hot and cold water control valve, high pressure cleaning nozzle with 2M hose, hot and cold water faucet.
- Automatic flushing table design.
- The height is adjusted by the electric mechanism, and the lifting stroke is 150mm.
- Waterproof leakage switch and socket
- Pipe and electrical line factory installation requires only a single point of connection at the project site.
- Size: 2500*800*750/900mm



SPINNE5



Product Name:

Autopsy Table.

Applications:

Anatomy, forensic examinations, medical colleges and hospitals, etc.

Product size:

L2180*W810*H900/1200mm; Special sizes can be customized.

Package:

Plywood case; packing size:2210*860*1030MM

Specification:

- Material: Full 304 stainless steel.
- Table hydraulic lift, height lift from 900 to 1200mm. It can be rotate 35 degrees.
- With ventilation system, and the adjustable air volume (need to connect with the user's ventilation system)
- Countertop/table concave design. At the bottom of the countertop, there are ventilation ducts and air volume control valves. The countertop is divided into two layers and the upper layer is two detachable perforated countertops. The lower floor is for ventilation and sewage discharge countertops.
- There is a sewage valve at the bottom of the lower deck. The sewage can flow into the water tank through the sewage valve.
- Removable writing device tray holder.

SPINNE4

Product Name:

Anatomical / autopsy table

Applications:

Mortuary; Anatomy, forensic examinations, medical colleges and hospitals, etc.

Product size:

Fixed model; L2400*W800*H750mm; Special sizes can be customized.

Lifting model; L2400*W800*H750 to 920mm; Special sizes can be customized.

Package:

Packed in plywood case; packing size:2430*790*970MM

Specification:

- Material: Full stainless steel 304.

Features:

- Circular arc connection within the table, no scaling, easy to clean.
- The anatomical table can be removed as a whole for easy cleaning.
- A corpse hair isolator is provided in the pool to prevent hair from entering the drain.
- There is a hand-operated faucet in the pool. The temperature of hot and cold water is adjustable (need to connect with the outside water heater)
- There is a mobile spray gun with adjustable water pressure that flushes the entire surface of the desk's anatomical table

Optional:

- Exhaust air design: Anatomy of the table top hole ventilation, (need to connect with the external exhaust fan) can be removed in time during the operation of harmful gases.
- There is a waterproof power socket on the side of the dissection table.
- The anatomical table is lifted up and down, using a new type of electric control lift, stable and reliable.
- Square base, fixed / lifting
- Size: length 2400* width 800* height (750/920) mm
- Power supply: AC 220V 50HZ



SPINNE5

Product Name:

Anatomical table / autopsy table

Applications:

Mortuary; Anatomy, forensic examinations, medical colleges and hospitals, etc.

Product size:

L1520*W760*H900mm; Special sizes can be customized.

Package:

Packed in plywood case; packing size:1550*810*1030MM

Specification:

- Material: Full stainless steel 304.
- Overall 304 stainless steel structure, surface grinding and wire drawing.
- All internal corner arc design for easy cleaning.
- 3 perforated panels for removable disassembly
- Ventilation: Downward suction.
- Aspirator and hot and cold water control valve, high pressure cleaning nozzle with 2M hose, hot and cold water faucet.
- Design for automatic table flushing
- The garbage handler can grind quickly
- Waterproof leakage switch and socket
- Pipe and electrical line factory installation requires only a single point of connection at the project site.
- Size: 1520*760*900mm

SPINNE7



Product Name:

Autopsy station

Applications:

Mortuary, Forensic examination, anatomy, medical colleges and hospitals, etc.

Product size:

L2130*W1060*H920mm; Special sizes can be customized.

Package:

Packed in plywood case; packing size: 2160*1110*1000MM

Specification:

- Material: Full 304 stainless steel.
- The frame is made of 30*30mm 304 stainless steel pipe.
- Whole welding, grinding drawing. With the wonderful making details.
- D=150mm mobile casters support brake.
- Countertop with a drain hole.
- The countertop is a concave design made of 304 stainless steel and is used with the stainless steel multi-function sink RD-118 or autopsy station.

SPINNEA (Working table thickness # 1.5mm) SPINNEB (Working table thickness # 2.0mm)

Product Name:

Small Dissecting Table

Applications:

Anatomy, forensic examinations, medical colleges and hospitals, etc.

Product size:

L1220*W760*H890mm; Special sizes can be customized.

Package:

Plywood case; Packing size: 1230*300*800MM

Specification:

- Material: Full 304 stainless steel.
- With anatomical support and writing board.
- With the Bottom drain and storage sewer.
- The Four-legged with wheels(with brakes) and for easy movement.
- The table with drain hole.



SPINNE-A1

Product Name:

Autopsy Assist table; Lifting type

Applications:

Anatomy, forensic examinations, medical colleges and hospitals, etc.

Product size:

L880*W580*H880/1040mm; Special sizes can be customized.

Package:

Plywood case; packing size: 910*1080*820MM



Specification:

- Material: Full 304 stainless steel.
- Adjustable height is from 880to 1040mm.
- 4 legs with the wheels/brakes
- With the stainless steel bucket.

THRIPSE



Product Name:

Autopsy Assist Equipment

Applications:

Laboratory, Forensic examination, anatomy, medical colleges and hospitals, etc

Product size:

L2900*W860*H1540mm; Special sizes can be customized

Package:

Packed in plywood case; packing size: 2930*890*1670MM

Specificati:

- Material: Full 304 stainless steele
- Large anatomical room, autopsy station supporting anatomical equipment.
- Function: It can be used with an anatomical cart. The dissection cart is moved to the sliding latch to connect with the sink. The operator can lock it from any side.
- The whole 304 stainless steel production, surface grinding and wire drawing processing.
- All internal corner arc design for easy cleaning. Pipes and electrical lines are installed in the factory, just point to point installation.
- Large drawers, LED lighting, garbage disposal and GFCI with waterproof cover socket.
- Luxurious stainless steel aspirator, power is water pressure vacuum, no power is more environmentally friendly and durable, water pressure control the size of the suction.
- Aspirator and hot and cold water control valve, high pressure cleaning nozzle with 3M hose, hot and cold water faucets.
- Ventilation duct design: Upper ventilation, or inner wall ventilation, needs to be connected to the user's ventilation system.
- Specifications: L2900W860*H1540mm

LAUS

Product Name:

Multifunctional washing table

Applications:

Laboratory, Forensic examination, anatomy, medical colleges and hospitals, etc.

Product size:

L810*W630*H1820mm; Special sizes can be customized

Package:

Packed in plywood case; packing size 1830*920*930MM

Specification:

- Material: Full 304 stainless steel.
- Multifunctional washing table.
- Made of high quality 304 stainless steel, surface grinding and wiredrawing.
- Configurations: faucets, showers, crushers and power strips.
- Luxurious stainless steel aspirator, powered by water pressure vacuum, no power is more environmentally friendly and durable water pressure to control the size of the suction



PHASMIDA1



Product Name:

Autopsy Table

Applications:

Mortuary; Anatomy, forensic examinations, medical colleges and hospitals, etc.

Product size:

L1500*W600*H850mm; Special sizes can be customized.

Package:

Packed in plywood case; packing size: 1550*680*980MM

Specification:

- Material: Full 304 stainless steel.
- Description: Animal Anatomy Table RD--1009
- The sinking anatomy bench is designed for autopsy and anatomy, surgery, and perfusion to provide a comfortable and hygienic platform.
- Standard design: There are 2 removable punching tops on the upper level for easy cleaning, with a sink and a hand-operated faucet or a sensor faucet. The lower deck is designed to sink more completely for overall sinking.
- Dissecting table specifications: L1500mm*W 600mm*H 850mm, the overall use of high-quality 304 stainless steel welding, all exposed welds have been seamlessly polished, the anatomical table is concave 1cm, to prevent liquid from flowing out of the table, the table is designed to move, you can free to move and other equipment supporting use.
- Sewage design: The lower end of the entire table top is inclined, connected to the sink drained to the sewer
- Cleaning design: High-pressure water gun cleaning, with 2 m anti-riot water pipe.
- Optional Center: 1. Anti-pollution vacuum cutting

PHASMIDA3

Product Name:

Autopsy table

Applications:

Mortuary; Anatomy, forensic examinations, medical colleges and hospitals, etc.

Product size:


L1980*W900*H850mm; Special sizes can be customized.

Package:

Packed in plywood case; packing size: 2010*950*980MM

Specification:

- Material: Full 304 stainless steel.
- Description: Convection-type animal dissection table
HZ-1003 Sinking dissection anatomy is designed to provide a comfortable and hygienic platform for autopsy and anatomy, surgery, and perfusion.
- The dissecting table functions to pull the air, odor, vapor, and aerosol from the person's face to a negative air source.
- Chemicals such as formaldehyde, ethanol, and xylenes can be controlled at lower than OSHA-allowed exposure limits, thereby providing a healthier and safer working environment.
- Standard design: There is one piece on the upper layer to remove the punching concave tabletop for easy cleaning, and the lower deck design for the overall sinking is more thorough.
- Dissector Bench Specifications: L1980*W900*H850mm The overall use of high-quality 304 stainless steel welding, all exposed welds are seamlessly polished, the anatomical table concave 1cm to prevent the liquid from flowing out of the table The table is designed with a double base, providing knees and legs Under the table, so that the user can choose the height of the stool seat or stand.
- Ventilation design: It can be pumped upwards through the ceiling connection, it can also be connected at the back through the wall, or downdraft through the ground to connect the outdoor exhaust fan, the exhaust pipe can be round or rectangular, the ventilation pipe has a manual air volume adjustment cutting You can freely adjust the air volume.
- Sewage design: the entire slope of the lower mesa is connected to the sewage outlet
- Cleaning design: high pressure water gun cleaning
- Power supply design: Equipped with power supply and leakage protection.

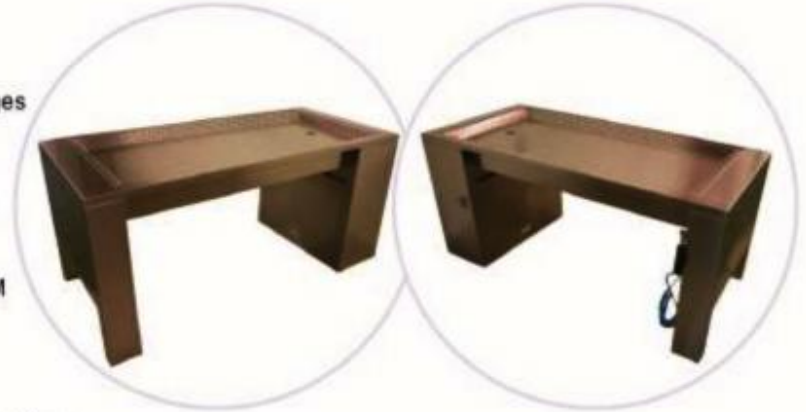

the sewer.

- Cleaning design: high pressure water gun cleaning
- Power supply design: equipped with power supply and leakage protection
- The dissecting table, which pulls down the air, odors, vapors, and aerosols, controls chemicals such as formaldehyde, ethanol, and xylene and can be kept below the allowable exposure limits of OSHA, thus providing a healthier and safer working environment
- Optional: Anti-pollution vacuum cutting