

AS-480

Auto Chemistry Analyzer

Technical specifications

General information

Machine type	Random access, open reagent & locked reagent (optional)
Test speed	Constant 400 tests/hour (Mono reagent or double reagent)
Test principle	Colorimetric method, turbidimetry
Test method	1 point end, 2 point end, Fixed-Time, Kinetics
Calibration type	Linear & nonlinear

Sample unit

Sample tray	120 sample positions
Sample tube	Micro cup & test tube & blood collection tube can be used. (Φ 12 -13)mm * (25 - 100)mm
Sample volume	2-30μL, step by 0.1μL
Sample probe	Liquid surface detection, timely tracking for liquid volume. Collision avoidance
Sample barcode	Sample barcode reader (optional)

Reagent unit

Reagent tray	90 reagent positions , Non-stop cooling system with peltier pad inside, 24 hours 2°C-14°C, open system, could accept any reagent.
Reagent volume	20-300μl, step by 1μl
Reagent needle	Liquid surface detection, timely tracking for liquid volume. Collision avoidance.
Reagent barcode	Reagent barcode reader (optional)

Reaction unit

Reaction tray	90 reaction cuvettes, new amorphous optical plastics reaction cups, perfect transmittance.
Reaction volume	150μl-330μl
Reaction time	10 minutes
Temperature	Peltier pad incubation system, reaction temperature (37°C ± 0.1 °C)

Mixing system

Mixer	paddle-type mixing, effective reducing of carry-over
-------	------------------------------------------------------

Optical system

Light source	Halogen lamp
Wavelength	340nm; 405nm; 450nm; 505nm; 546nm; 578nm; 630nm; 700nm(4 more options)
Absorbance range	0-4.0 Abs
Resolution	0.0001 Abs

Operating system

Operation system	Windows XP, Windows 7, Windows 8, Windows 10
Data processing	Automatic calibration, test combination, effective time management of reagents, serum index, whole process detection, linear expansion of enzyme, blank deduction, dirty cuvette memory , cross infection prevention, patient information memory and association input, automatic report audit, data fuzzy query, report statistics and printing, reference range grading, alarm information classification, user operation privilege classification management.
Report printing	6 formats optional, support custom pattern
interface	TCP/IP network interface

Working conditions

Power supply	AC 220V, 50Hz; AC 110V, 50/60Hz, 1000VA
Water consumption	15L/hour(maximum)
Dimension	1100mm * 765mm * 1150mm(L*W*H)
Weight	165Kg



AS-480

Auto Chemistry Analyzer

400 Tests/Hour, Double Reagent

BIOELAB

THE
DRAKE
BUSINESS GROUP
CORPORATION



Drake Marketing and
Equipment Corporation



Drake Calibration
Laboratory Philippines



Erpel Konstrukt
Konzern Corporation



Urossorsa Marketing and
Equipment Corporation

(02) 8463 4943 | (02) 8463 5120
(+63) 939 911 9277

1052 Benigno Aquino Avenue, Sitio 2
Bahay Pare, Candaba, Pampanga
Philippines, 2013

176 Benigno Aquino Avenue, Sitio 2
Bahay Pare, Candaba, Pampanga
Philippines, 2013

information@drakephilippines.com

AS-480


Auto Chemistry Analyzer

- Constant 400 test / hour with double reagent
- 90 reagent positions
- Automatic liquid level sensing, anti-collision protection.
- Integrated full closed optical system, multi wavelength simultaneous detection.
- New amorphous optical plastics reaction cups, perfect transmittance.
- Maintenance free semiconductor solid state refrigeration
- Friendly surface and easy to operation
- 6-stage, 2 times automatic washing system , Alkaline cleaning solution



Multifunctional precision sampling probe

- External & internal mirror polish, external & internal probe washing
- Dedicated sampling probe equipped with sensitive liquid sensor, timely feedback of reagent & sample residuals
- Collision protection, automatic probe depth adjusting

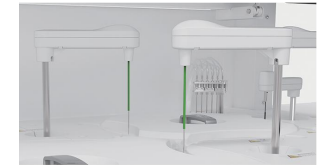

Excellent mixer design

- Teflon coating mixer. No water dropping (reducing carry-over)
- Excellent mixing effect with standard mixing procedure



High precision ceramic syringe

- Permanent ceramic piston
- Accurate dispensing as low as 0.1 μ l



Reagent & sample tray

- 90 reagent positions
- 120 sample positions, micro cup & test tube can be used
- Non-stop cooling system with peltier pad inside, 24 hours 2 $^{\circ}$ C-14 $^{\circ}$ C
- Bar code reader(optional)



Reaction tray

- 90 reaction cuvettes
- Reaction volume as low as 150 μ l
- Stable & accurate temperature(37 \pm 0.1 $^{\circ}$ C) for reaction



Advanced operating system

- User friendly software interface simple, easy-to-operation
- Sample automatic dilution(decrease, increase, normal)
- Reaction result automatic checking and warning, calibration result automatic checking, warning alarm log
- Linear & nonlinear calibration,Multi-points calibration up to 8 points
- Quality control apply to westgard, L-T plot, cumulative

