


## HEMOGLOBIN ANALYZER

- 7 $\mu$ L Blood
- HB/HCT
- $\leq$  3s testing time



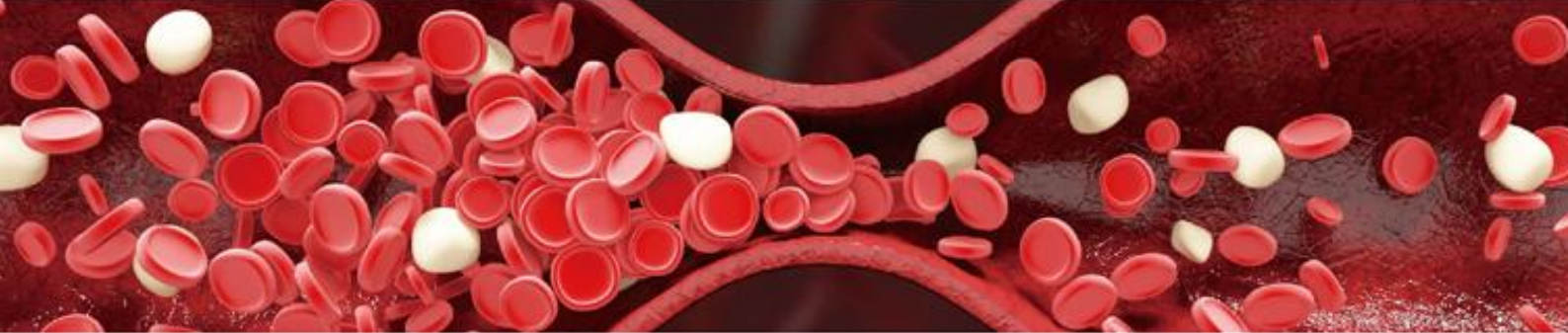


Parameters: Hb and Hct ●

Measurement range: ●  
0–25.6 g/dL

Shelf life: 24 months sealed, ●  
3 months unsealed

Packing: 50 pcs/bottle ●

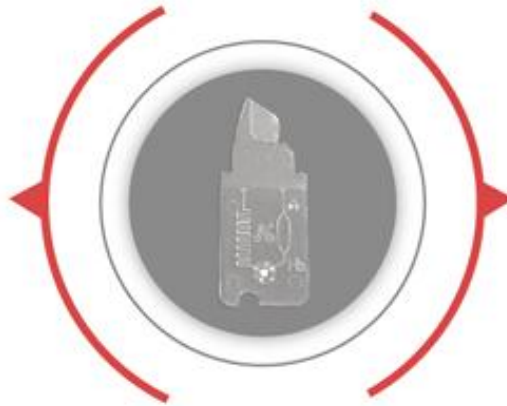


● Reagent free

● Results: ≤3s

● Blood sample: 7μL,  
Capillary venous whole blood

● Storage Condition:  
Room temperature

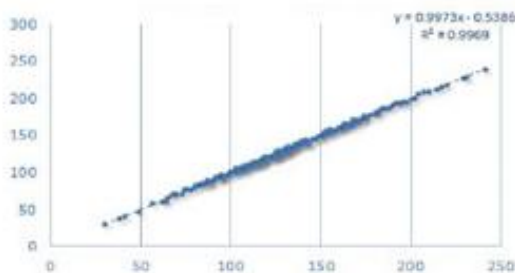

There is sufficient evidence to conclude that there is a significant linear relationship between the test results from Konsung H7 and laboratory analyzers.

## Analyzer Specifications

Principle	Microfluidic and spectrophotometry with scattering compensation technology	Operation environment	5°C ~ 35°C ; ≤85% RH
Calibration	Factory calibrated needs no further calibration	Interface	USB, Bluetooth, wifi, printer
Power supply	AC100 ~ 240V, 20VA DC 5V, 1A Built-in lithium battery	Measurement	130mm × 82mm × 31.5mm
Data storage	2000 test results	Weight	500g
Accuracy	CV ≤ 1.5%	Optional	Bluetooth printer APP

THE  
**DRAKE**  
BUSINESS GROUP  
CORPORATION

(02) 8463 4943 | (02) 8463 5120  
(+63) 939 911 9277



Drake Marketing and  
Equipment Corporation



Drake Calibration  
Laboratory Philippines



Erpel Konstrukt  
Konzern Corporation



Urossorsa Marketing and  
Equipment Corporation

1052 Benigno Aquino Avenue, Sitio 2  
Bahay Pare, Candaba, Pampanga  
Philippines, 2013

176 Benigno Aquino Avenue, Sitio 2  
Bahay Pare, Candaba, Pampanga  
Philippines, 2013

information@drakephilippines.com